


GE iDXA enCORE v. 17 & 18

Total Body Densitometry – save as *ID_visit#_Bone*

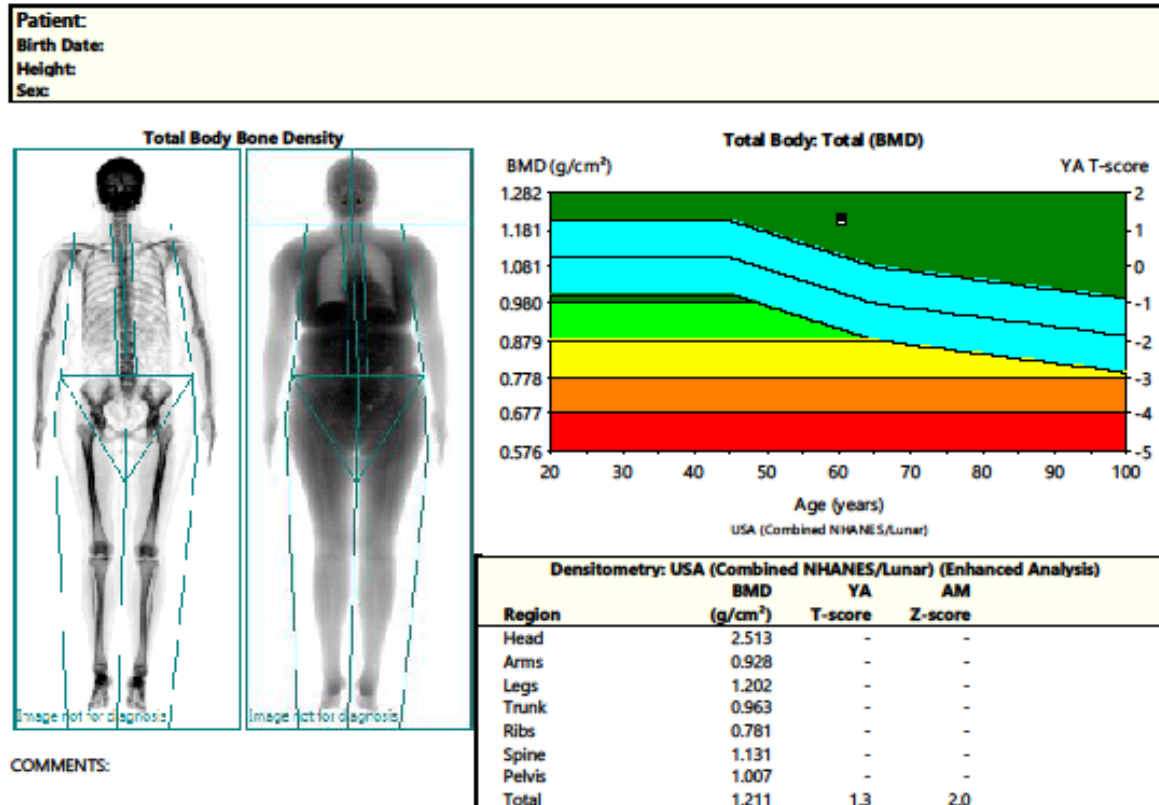

Statistically 68% of repeat scans fall within 1SD (± 0.010 g/cm³ for Total Body Total); USA (Combined NHANES (ages 20-30) / Lunar (ages 20-40)) Total Body, Female Reference Population (v113); Matched for Age, Sex, Weight (females 25-100 kg), Ethnic:

Date created: 02/01/2020 9:02:27 AM 17 [SP 4]; Filename: 14ej4q4jy.mab; Total Body: 100;0.19;153.85;15.6 0.00-1.00 2.40;3.04 13.2;5%Fat=49.2%; 0.00;0.00 0.00;0.00; Scan Mode: Standard; 3.0 μ Gy




GE Healthcare

Page: 1 of 1

Lunar iDXA
ME+212290

GE iDXA enCORE v. 17 & 18

Total Body Composition – save as *ID_visit#_BC*

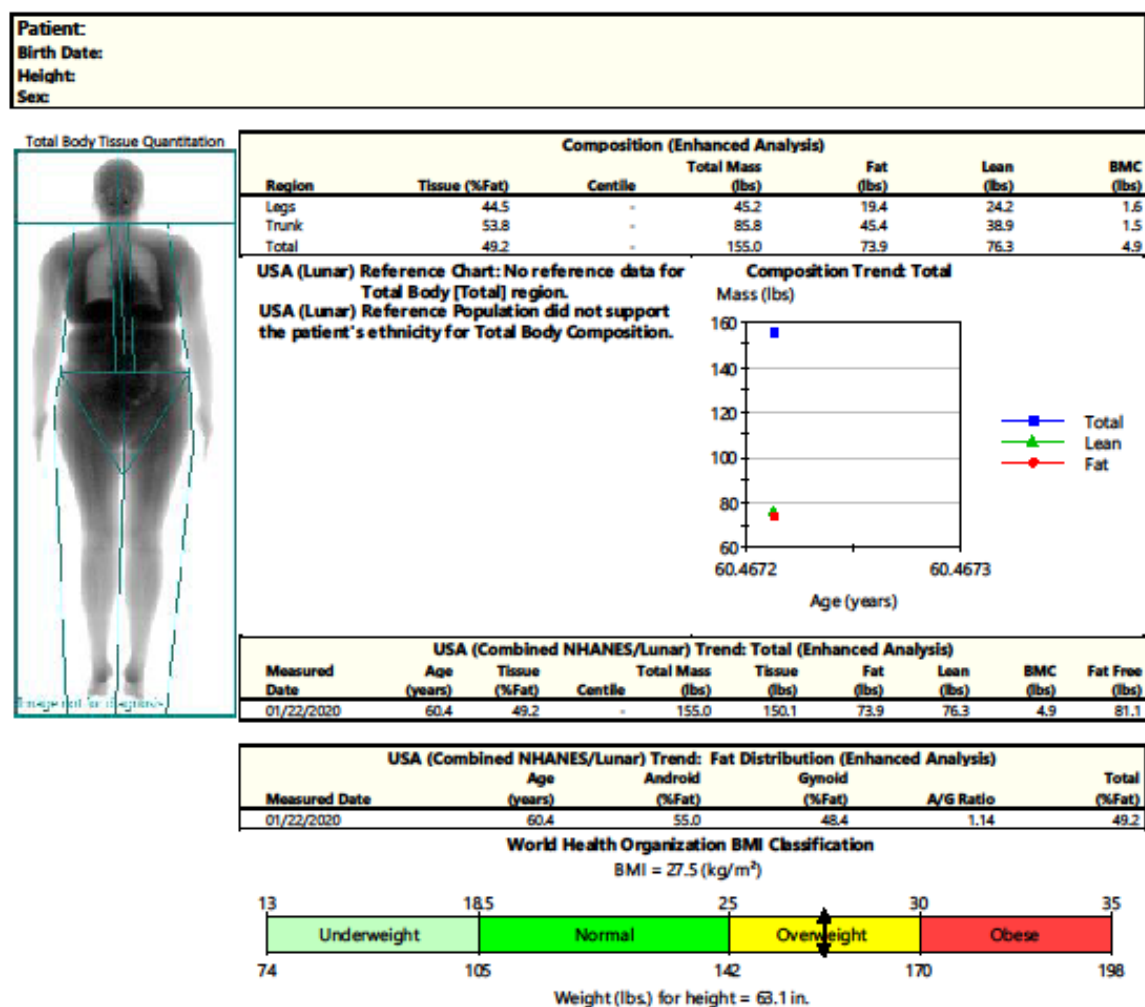

COMMENTS:

Statistically 68% of repeat scans fall within 1SD ($\pm 0.4\%$ Fat, ± 0.33 lbs. Tissue Mass, ± 0.62 lbs. Fat Mass, ± 0.68 lbs. Lean Mass for Total Body Total); USA (Lunar) Total Body Composition, Female Reference Population (v113); Composition Matched for Age, Sex
Date created: 02/01/2020 9:52:28 AM 17 [SP 4]; Filename: 14ej4q4pyzmb; Total Body: 100,0.19:153.85:15.6 0.00:-1.00 2.40:3.04 13.2%Fat=49.2%; 0.00:0.00 0.00:0.00; Scan Mode: Standard; 3.0 uGy



GE Healthcare

Page: 1 of 1

Lunar iDXA
ME+212290

GE iDXA enCORE v. 17 & 18

Total Body Composition Ancillary – save as *ID_visit#_BCExt*

Patient:	Referring Physician: (not specified)
Birth Date:	Patient ID:
Height:	Measured:
Sex:	Analyzed:

BODY COMPOSITION: Total Body (Enhanced Analysis)

Region	Tissue (%Fat)	Region (%Fat)	Tissue (lbs)	Fat (lbs)	Lean (lbs)	BMC (lbs)	Total Mass (lbs)
Arms	0.0	0.0	0.0	0.0	0.0	0.0	0.0
Legs	60.8	58.8	9.2	5.6	3.6	0.3	9.5
Trunk	62.8	60.6	14.9	9.4	5.5	0.6	15.5
Android	63.3	61.0	1.7	1.1	0.6	0.1	1.8
Gynoid	64.0	61.6	3.4	2.2	1.2	0.1	3.5
Total	59.9	57.9	26.0	15.5	10.4	0.9	26.9

Fat Mass Ratios:

Trunk Fat Mass/Total Fat Mass	Legs Fat Mass/Total Fat Mass	Limbs Fat Mass/Trunk Fat Mass
0.60	0.36	0.60

Estimated Visceral Adipose Tissue

Volume	Mass	Area
25.22 in ³	0.86 lbs	28.50 in ²

Estimated Subcutaneous Adipose Tissue


Volume	Mass	Area
11.00 in ³	0.38 lbs	12.44 in ²

Statistically 98% of repeat scans fall within 1SD ($\pm 0.4\%$ Fat, ± 0.33 lbs. Tissue Mass, ± 0.62 lbs. Fat Mass, ± 0.68 lbs. Lean Mass for Total Body Total)

Date created: 10/20/2020 12:35:05 PM 18 [SP 2] Filename: qpgliq4jeo.melt; Total Body: 100/0.19/153.85/15.6 0.00/-1.00 2.40/3.04 13.0/%Fat=59.9%; 0.00/0.00 0.00/0.00; Scan Mode: Standard; 3.0 μ Gy; 1.20 cGy/cm²

GE iDXA enCORE v. 17 & 18

Lumbar Spine – save as *ID_visit#_LS*

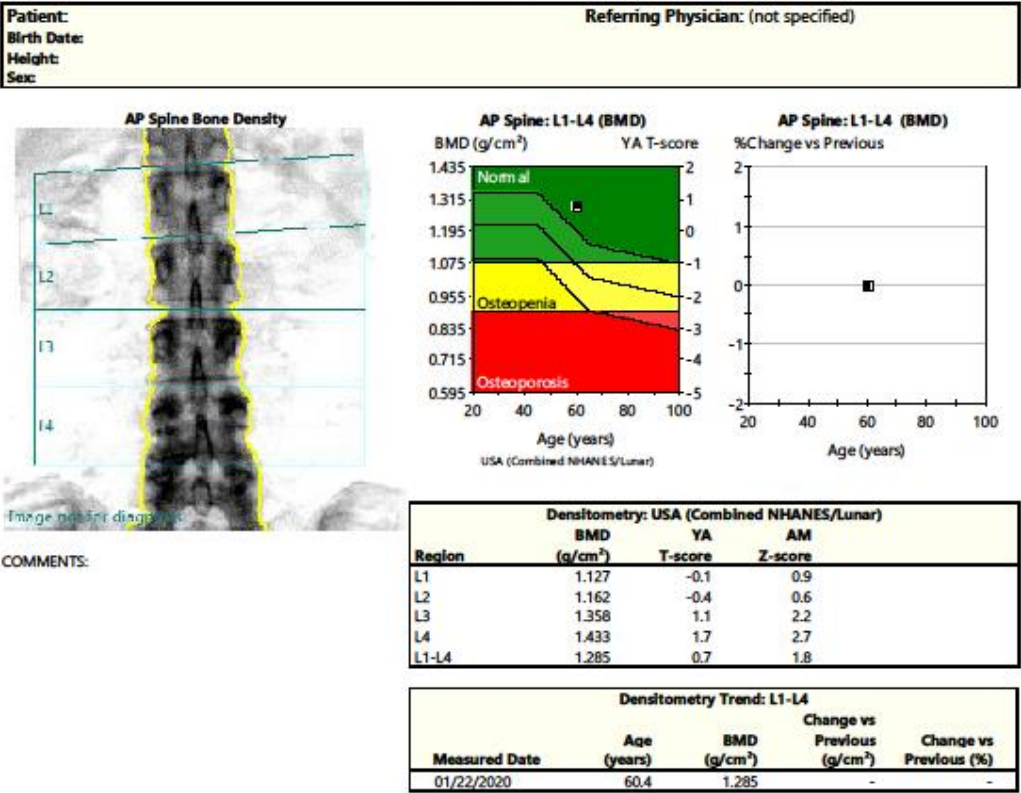

(*) Indicates significant change based on 95% confidence interval. (LSC= 0.028 g/cm³ for AP Spine L1-L4); Statistically 68% of repeat scans fall within 1SD (± 0.010 g/cm³ for AP Spine L1-L4); USA (Combined NHANES (ages 20-30) / Lunar (ages 20-40)) AP Spine, Female Reference Population (v113); Matched for Age, Sex, Weight (females 25-100 kg), Ethnic; World Health Organization - Definition of Osteoporosis and Osteopenia for Caucasian Women: Normal = T-score at or above -1.0 SD; Osteopenia = T-score between -1.0 and -2.5 SD; Osteoporosis = T-score at or below -2.5 SD; (WHO definitions only apply when a young healthy Caucasian Women reference database is used to determine T-scores)


Date created: 02/01/2020 9:59:41 AM 17 (SP 4); Filename: 1faj4qjpy.mec; AP Spine; 100.25050.0050.0 0.00:10.74 0.30x0.25 23.5;Hfat=41.7%; 0.00:0.00 0.00:0.00; Scan Mode: Standard;OneScan: 146.0

µg

GE Healthcare

GE iDXA enCORE v. 17 & 18

Dual Hip – save as *ID_visit#_DH*

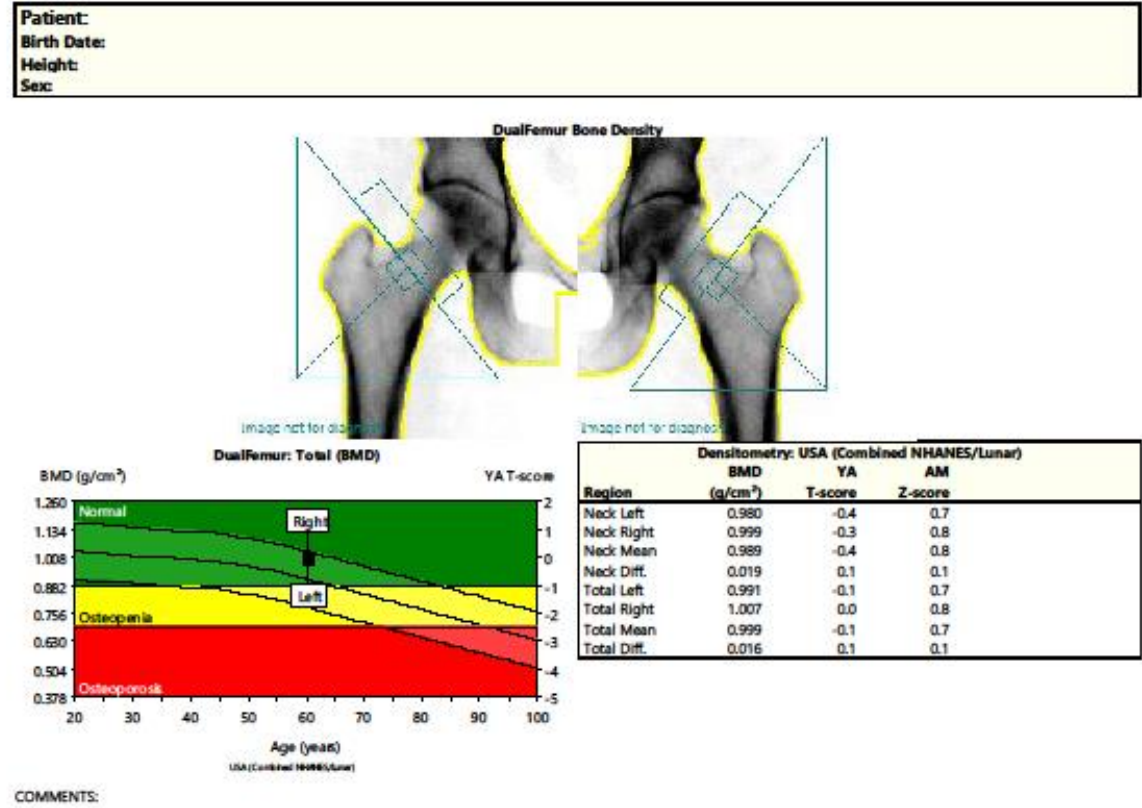

GE iDXA enCORE v. 17 & 18

R or L Radius – save as *ID_visit#_Rad*

** if you scan both radii you will only get data from a single scan, place the other side scans in a separate input folder and repeat the process and simply cut and past the data into your main data document.*

